



# Mobile Holder for Image Receptors

Cat# MCH-33C

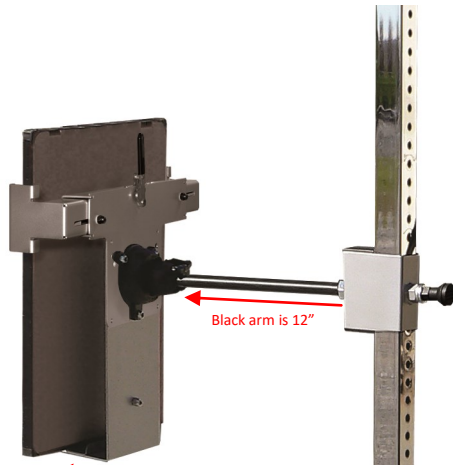

Front of cassette is  
18" from pole


Vertical Range from bottom of  
image receptor: 1" off Floor to 54"  
allowing ample room for imaging  
at several height levels.

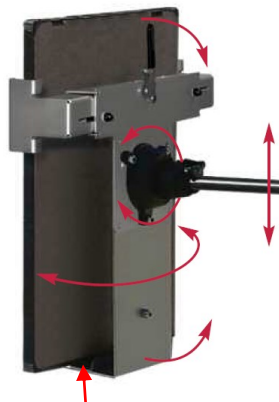

1.5" Deep Channel Accepts Grids and Caps  
and will accommodate most receptor sizes  
like: Agfa, AllPro, Carestream, Canon, Dell,  
Fuji, GE, Kodak, Konica, Quantum, Medlink,  
Rayence, Samsong, Shimadzu, Toshiba,  
Varian, and Vieworks

## Fully Adjustable Tilt & Rotation Abilities:

- \* Straight Vertical +/- 90 degrees
- \* Incline Forward / Backward
- \* Horizontal rotation +/- 180 degrees
- \* Full Rotate or Oblique Tilt
- \* Face Up
- \* Face Down
- \* 360 degrees - Anywhere in Between

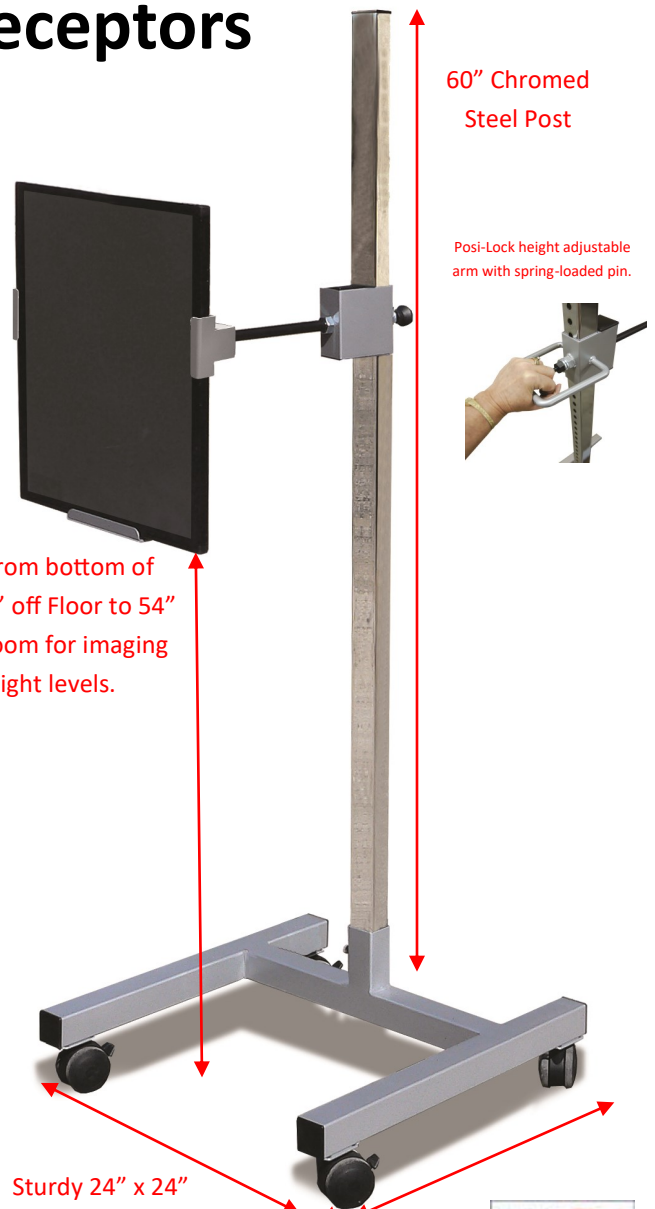

60" Chromed  
Steel Post

Posi-Lock height adjustable  
arm with spring-loaded pin.

Sturdy 24" x 24"  
Broad H-base  
prevents tipping.

Heavy Duty locking  
3" casters for easy  
manipulation.



Image Receptor Obliqued at Angle

Image Receptor Face Up

Image Receptor Tilted Backwards

Radiographs can be taken from any angle. Tray can be positioned face up, face down, and anywhere in between to maximize proper receptor positioning. Perfect for ER, OR, and x-ray procedures such as myelograms, decubitus studies, and three way acute abdomen exams.

- \* Bilateral Clamps securely holds receptors in place. Holds up to 20.5" wide flat imaging receptors with No Height Limitation.
- \* Will hold receptors up to 12lbs.
- \* Easy to clean in a "OR" or "Sterile" setting.
- \* Overall equipment weight is approx 54 lbs.
- \* Manufactured in Nashville, Tennessee USA
- \* 14 Day Lead Time
- \* One Year Warranty
- \* Ships in 3 cartons: Pole 12lbs, base 23 lbs , & arm 19 lbs.