

Evolve Casework and Cabinets

Storage Solutions for Healthcare



More than ordinary casework.

Evolve cabinets and casework make the work of storing and accessing supplies efficient and organized, and they do so with a clean design that offers enhanced aesthetics and construction details that support infection prevention measures. The ROI of Evolve modular casework translates to depreciation benefits, faster installation, and adaptability to change.

The Evolve line includes preconfigured cabinets, cabinets with FlexCell and a choice of interior components, upper and lower cabinets, catheter cabinets, and scope cabinets offered in several sizes and a range of material and finish options, including our proprietary AireCore construction, which is resistant to moisture, bacteria, and corrosion.



Products shown, Left to right: Cabinet with accessory pack (SEB2736TGAP), Cabinet with divided shelves (SEB2736TGDS), Cabinet with FlexCell and split center column (SEB2736TGSC). Accessories sold separately



A casework approach that creates a cohesive line of modular cabinetry.

Flexibility in function and style. Evolve casework and its FlexCell system of shelves, dividers, trays, and baskets let you arrange the interior configuration of your cabinets for specific types and sizes of supplies and easily reconfigure as needs and supplies change. Cabinets with FlexCell can be specified with glass or solid hinged doors and roll-top doors.



Finish the sides of cabinets with optional finished ends that match the cabinet's surface material or add an infill panel to the end of a cabinet run to fill gaps.




Sloped tops can be specified on new cabinets or retrofit after purchase.

Products shown, Left to right: Cabinet with accessory pack (SEC2736TGAP) and sloped top (SEC2736THST), Upper cabinets with adjustable shelf (SEC1336UG) and sloped top (SEC1336UHST), Base cabinet with FlexCell (SEC2736BSNC) and countertop (SEC2736CT), Catheter cabinet (SEC2736TGCC) with sloped top (SEC2736THST) in top image and divided shelves (SEC2736TGDS) in bottom image

Eliminate wasted space. An Evolve cabinet with FlexCell interior and a variety of trays and baskets has a storage capacity of 31.5 linear feet. Compare that to the 14.4 linear feet of storage capacity in a cabinet with adjustable shelves.*




*Comparison of two 36"w, 27"d, and 84"h cabinets with hinged doors. The cabinet with FlexCell has a center column and 60 cells (30 in each compartment). The other cabinet has four full-width adjustable shelves.

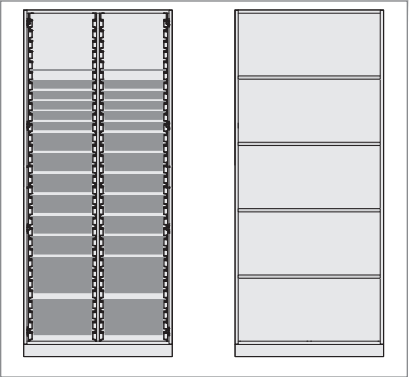

LED lights in cabinets provide visibility to supplies in hybrid ORs.

Products shown, Top, Left to right: Cabinet with accessory pack (SEC2736TRTAP), Cabinet with divided shelves (SEC2736TRTDS), Cabinet with FlexCell and split center column (SEC2736TRTSC). Accessories sold separately. Bottom, Left to right: Cabinet with accessory pack (SEB2736TGAP), Cabinet with divided shelves (SEB2736TGDS), Cabinet with split center column (SEB2736TGSC). Accessories sold separately

The details make the difference. Evolve casework is designed and engineered for healthcare environments.

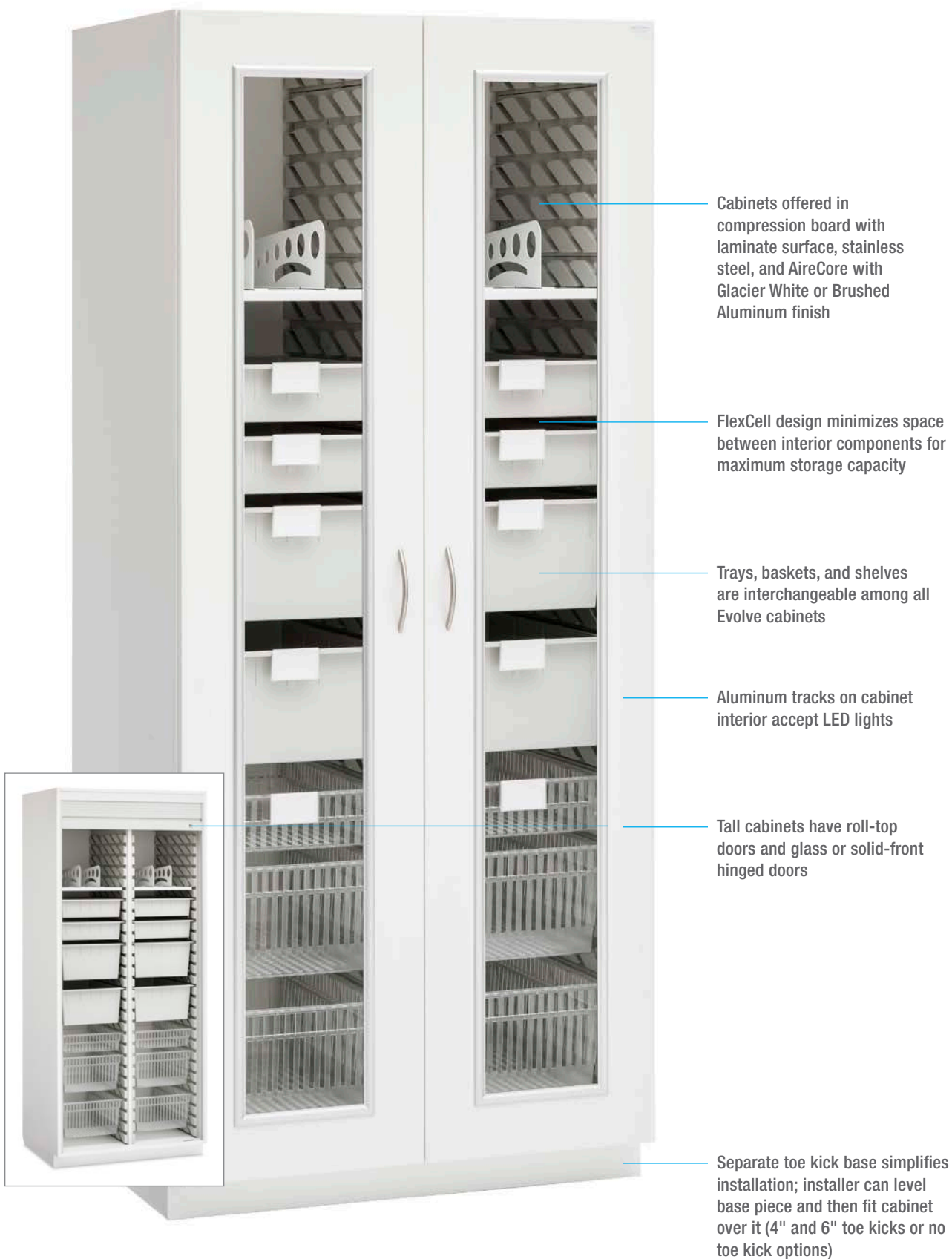

Products shown: Cabinet with accessory pack and glass door (SEA2736TGAP) in large image and roll-top door (SEA2736TRTAP) in smaller image



Upper cabinets have glass and solid-front hinged doors

Seamless edge banding eliminates glued and cracked edges



Supporting infection prevention initiatives. Our proprietary AireCore construction wraps an impervious material around an interior structure that is durable and lightweight. Combine AireCore with our seamless edge banding and you have a product that is resistant to moisture, bacterial growth, and corrosion. AireCore offers design choices, from roll-top or hinged doors to surface materials and customization.



Countertops for base cabinets are compression wood and laminate or stainless steel

Self-closing door with concealed hinges

Base cabinets have roll-top doors and solid-front hinged doors

Products shown: Upper cabinet with adjustable shelf (SEA1336UG), Base cabinet with FlexCell (SEA2736BSC) and countertop (SEC27CT). Accessories sold separately

Specialized Evolve cabinets are designed for the equipment and requirements of cardiac catheter, interventional radiology, and endoscopy departments, among a hospital's busiest.

Catheter cabinets, installed individually or in a run, efficiently and safely store catheters and related equipment and supplies. Scale the size of your cabinets to the square footage of the department to hold more catheters for use in shared spaces or fewer in individual procedure rooms.



Each inline hook holds 10-15 catheter packages and can be adjusted horizontally

Self-closing glass and solid-front doors feature concealed hinges

Aluminum channel offers durable, finished detailing on door

AireCore with Brushed Aluminum finish is a smart alternative to stainless steel storage in surgical suites and procedure rooms

Products shown: Catheter cabinet (SEB1936TGCC)



Boxed catheter cabinets with pull-out boxes and dividers offer easy access and organize sizes and style of catheters.



Catheter storage products, designed with FlexCell interiors, maximize storage space and keep supplies accessible and easily identified.

Products shown: Top, Left to right: Boxed catheter cabinet with glass door (SEC1936TGBC) and with roll-top door (SEC1936TRTBC). Bottom, Left to right: Boxed catheter cabinet with left-hinged glass door (SEC2719TGLBC), Catheter cabinet (SEC2736TGCC), Cabinet with FlexCell and center column (SEC2736TGC), Cabinet with divided shelves with right-hinged glass door (SEC2719TGRDS), all with sloped tops. Accessories sold separately

Evolve scope cabinets provide ventilated, clean, cushioned, and dust-free environments for endoscopes. Scope cabinets come preconfigured and in sizes that accommodate different scope lengths and quantities. Scope managers and cord managers keep contents secure and organized.



Vents circulate air



Optional HEPA filtered fan unit provides positive air flow

Padded rear wall protects scope optics



Adjustable cord managers keep endoscopes organized, secure, and elevated



Removable drip tray eases cleaning



Products shown: Scope cabinet for 5 scopes with roll-top door (SEA1927SRTS) in small image and for 8 scopes with glass doors (SES1936SGS) in large image



Cabinets can be laminated to match surrounding décor.



Cabinet stores 5 to 18 scopes up to 83" in length for cabinets with hinged doors and 76" for cabinets with roll-top doors.

Adjustable single and double scope managers adapt to different lengths of scopes.



The Evolve scope drying cabinet includes all the features of the scope cabinet but also offers the drying unit that provides automated, consistent, and continuous delivery of HEPA-filtered, pressurized air into all endoscope channels to keep scopes dry and avoid reprocessing. The easy-to-operate system features an automated, 10-minute drying cycle per scope and continuous venting cycle for overnight storage. It has self-adjusting, high-pressure regulation, automatic shut-off, and a 99.9% retention rate of particles less than 0.3µm.

The drying unit can be ordered separately and field-installed in a new or previously purchased Evolve scope cabinet with AireCore or stainless steel construction.

Products shown, Left to right: Scope cabinet for 8 scopes (SEC1936SGS), Upper cabinet with adjustable shelf (SEC1336UG), Base cabinet with FlexCell (SEC2736BSNC) and countertop (SEC27CT), Scope drying cabinet (SEA1927SGRSF)

Configure Evolve cabinets for high-density storage without taking up valuable floor space.



Preconfigured options include cabinets with divided shelves.

FlexCell is fastener free—no tools are required to reconfigure interiors.

2"-, 4"-, and 8"-high trays match the size of supplies to the size of the tray and reduce wasted vertical space.



Stackable trays make transporting restocked trays efficient.



Products shown, Left to right: Cabinet with divided shelves (SEC1936TGDS), Cabinet with FlexCell and offset column (SEC1940TGRC), 2" Trays (ST2). Accessories sold separately



Trays and baskets offer moderate and gravity feed tilts. Dividers are adjustable and removable so you can organize and reorganize tray contents.



Clear baskets help to identify contents quickly; perforated bottoms prevent dust build-up.



Manage access to high-cost supplies using proximity locks. An added level of security is available with dual protection via proximity lock and keypad code. Proximity locks also give staff quick access to supplies, especially important in emergency departments.



Quickly locate and restock supplies with the help of clip-on labels.

Products shown, Top, Left to right: 2" Trays with dividers (STD2), 4" Basket (SB4). Bottom, Left to right: Keypad lock (SLE), Keypad with proximity lock (SLD), Label holders (SLH35)


Accessory packs offer a streamlined way to outfit cabinet interiors of any size.




Products shown, Left to right: Cabinet with accessory pack, 36" w x 27.25" d (SEC2736TRTAP) and 26.75" w x 19.25" d (SEC1927TRTAP)




2" Tray with Dividers (STD2)

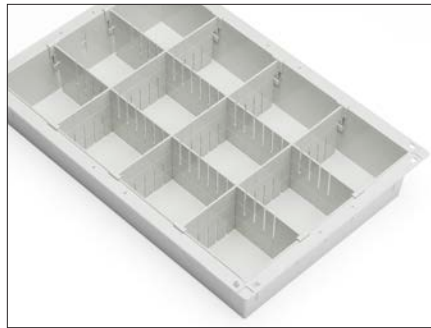

4" Tray with Dividers (STD4)

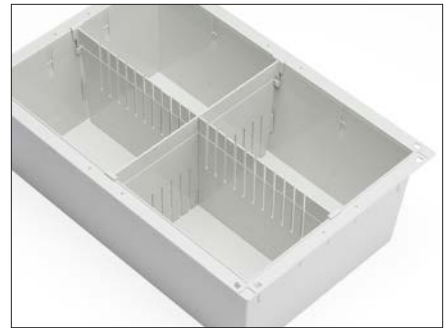


8" Tray with Dividers (STD8)

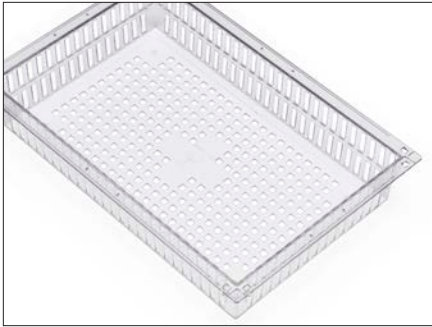

4" Basket (SB4)



8" Basket (SB8)

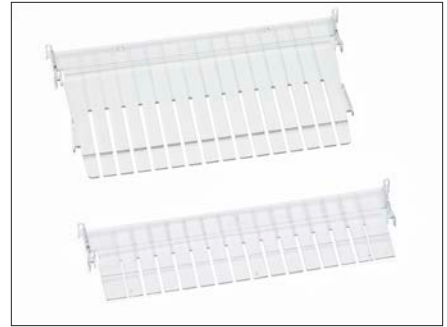

4" and 8" Short Basket Dividers (SBD4, SBD8)



2" Clear Tray with Dividers (SCT2)



Label Holders (SLH35)



Metal Shelf (S27S1, S19S1)



Multipurpose Shelf (S27M1, S19M1)
with 3-Cath Manager (S9CM, S7CM)



Lockable Drawer for single compartment (S27D1, S19D1)



Multipurpose Shelf (S27M1, S19M1)
with 2 Dividers (S27D1, S19D1)



Multipurpose Shelf (S27M2, S19M2)
with 7-Cath Manager (S9CM, S7CM)



Lockable Drawer for double compartment
(S27D2, S19D2)



Multipurpose Shelf (S27M2, S19M2)
with 5 Dividers (S27D1, S19D1)

Evolve cabinets and casework create a flexible storage system that organizes supplies and equipment throughout a hospital.

Evolve cabinets are offered in four widths and three depths; each cabinet type has one height. Scope cabinets are offered in additional heights to accommodate varying scope lengths. Preconfigured divided shelf, accessory pack, catheter, and boxed catheter cabinets are available. Cabinet options include sloped tops, keyed and keyless locks, toe kicks or no toe kick, and interior accessories.



Evolve wall-mounted
scope cabinet
26.75"w x 48"h

Evolve scope cabinet
26.75"w x 60"h

Evolve scope cabinet
26.75"w x 92"h

Evolve scope cabinet
36"w x 92"h

Evolve scope cabinet
40"w x 92"h

Evolve scope cabinets

Evolve wall-mounted cabinet

The 26.75"-wide, 19.25"-deep, 48"-high wall-mounted cabinet includes a lockable counterbalanced roll-top door or lockable tempered glass door that can be hinged on the left or the right. The cabinet is available with a scope or ultrasound cabinet interior. The ultrasound cabinet is available with glass door only.

Exterior dimensions:

26.75"w x 19.25"d x 48"h

Interior configuration:

Scope cabinet holds 5 scopes; scope hanging length is 40" maximum (hinged door) or 35" maximum (roll-top)
Ultrasound cabinet holds 5 probes

Evolve 26.75"w x 60"h scope cabinet

The 26.75"-wide, 19.25"-deep, 60"-high scope cabinet includes a lockable counterbalanced roll-top door or lockable tempered glass door that can be hinged on the left or the right. Cabinet has a 4"-high toe kick.

Exterior dimensions:

26.75"w x 19.25"d x 60"h

Cabinet with 6" toe kick is 62"h

Interior configuration:

Holds 5 scopes; scope hanging length is 51" maximum (hinged door) or 44" maximum (roll-top)

Evolve 26.75"w x 92"h cabinet

The 26.75"-wide, 19.25"-deep, 92"-high cabinet includes a lockable counterbalanced roll-top door or lockable solid-front or tempered glass door that can be hinged on the left or the right. The cabinet is available with a scope or TEE probe cabinet interior. Cabinet has a 4"-high toe kick.

Exterior dimensions:

26.75"w x 19.25"d x 92"h

Cabinet with 6" toe kick is 94"h

Interior configuration:

Scope cabinet holds 5 or 10 scopes; scope hanging length is 83" maximum (hinged door) or 76" maximum (roll-top)
TEE probe cabinet holds 3 probes

Evolve 36"w x 92"h cabinet

The 36"-wide, 19.25"-deep, 92"-high cabinet includes a lockable counterbalanced roll-top door or lockable solid-front or tempered glass doors. The cabinet is available with a scope or TEE probe cabinet interior. Cabinet has a 4"-high toe kick.

Exterior dimensions:

36"w x 19.25"d x 92"h

Cabinet with 6" toe kick is 94"h

Interior configuration:

Scope cabinet holds 8 or 16 scopes; scope hanging length is 83" maximum (hinged door) or 76" maximum (roll-top)
TEE probe cabinet holds 5 probes

Evolve 40"w x 92"h cabinet

The 40"-wide, 19.25"-deep, 92"-high cabinet includes solid-front or tempered glass doors. Cabinet has a 4"-high toe kick.

Exterior dimensions:

40"w x 19.25"d x 92"h

Cabinet with 6" toe kick is 94"h

Interior configuration:

Holds 9 or 18 scopes; scope hanging length is 83" maximum (hinged door)

Product numbering system for Evolve scope cabinets



Product Shown:
Scope cabinet for 8 scopes
SES1936SGS

S	E	S	1	9	3	6	S	G			S
Product Line	Material	Depth	Width	Cabinet Type	Door Type	Hinge Location	Cabinet Configuration				

Product Line
SE – Evolve
Material
A – AireCore
B – AireCore with Brushed Aluminum
C – Compression Board
S – Stainless Steel
Depth (rounded to 0.25")
19 – 19"
Width (rounded to 0.25")
27 – 27"
36 – 36"
40 – 40"

Cabinet Type
S – Scope
Door Type
G – Glass
RT – Roll-top
S – Solid-front
Hinge Location
R – Right
L – Left

Cabinet Configuration
S – Single Scope Holder, for 92"h Scope Cabinet
D – Double Scope Holder, for 92"h Scope Cabinet
W – Wall-mounted Scope Cabinet
B – 60"h Bronchoscope Cabinet

Products shown: Scope cabinet for 8 scopes (SES1936SGS)

Evolve's flexibility extends to material choices.

AireCore construction with thermally fused, seamless edge banding is resistant to water, chemicals, and bacteria. AireCore with a Brushed Aluminum finish offers the look of stainless steel and the same nonporous surface but without fingerprints and smudges, making it a low-maintenance option. AireCore can also be finished in Glacier White and standard or custom laminate. Stainless steel and compressed board with laminate are also available.



Evolve upper and base cabinets

Base cabinet

The 34.5"-high base cabinet comes in 3 widths and 2 depths. Cabinet has a lockable solid-front or roll-top door and 4"-high toe kick. Order countertop separately.

Exterior dimensions:

- 18.75"w x 27.25"d x 34.5"h
- 26.75"w x 19.25"d x 34.5"h
- 36"w x 27.25"d x 34.5"h
- Cabinet with 6"h toe kick is 36.5"h; cabinet with no toe kick is 30.5"h

Interior configuration:

- 11 FlexCells, 22 in cabinet with center column (hinged door)
- 8 FlexCells, 16 in cabinet with center column (roll-top door)
- 36"-wide cabinet available with or without center column

Upper cabinet

This 13"-deep, 30"-high upper cabinet comes in 3 widths and has 1 adjustable full-width shelf. Door options include tempered glass or solid-front doors.

Exterior dimensions:

- 18.75"w x 13"d x 30"h
- 26.75"w x 13"d x 30"h
- 36"w x 13"d x 30"h

Interior configuration:

- Shelf adjusts in 1.25" increments

Evolve 19.25"d cabinets

26.75"w cabinet

The 26.75"-wide, 19.25"-deep, 84"-high includes a counterbalanced roll-top door or solid-front or tempered glass door that can be hinged on the left or the right. Cabinet has a 4"-high toe kick.

Exterior dimensions:

- 26.75"w x 19.25"d x 84"h
- Cabinet with 6" toe kick is 86"h; cabinet with no toe kick is 80"h

Interior configuration:

- 30 FlexCells (hinged door)
- 27 FlexCells (roll-top door)

36"w cabinet

The 36"-wide, 84"-high cabinet is available with no center column. Door choices include a counterbalanced roll-top door or solid-front or tempered glass doors. Cabinet has a 4"-high toe kick.

Exterior dimensions:

- 36"w x 19.25"d x 84"h
- Cabinet with 6" toe kick is 86"h; cabinet with no toe kick is 80"h

Interior configuration:

- 30 FlexCells (hinged door)
- 27 FlexCells (roll-top door)

40"w cabinet

The 40"-wide, 19.25"-deep, 84"-high cabinet includes a wide and narrow compartment. Door choices include solid-front or tempered glass doors. Cabinet has a 4"-high toe kick.

Exterior dimensions:

- 40"w x 19.25"d x 84"h
- Cabinet with 6" toe kick is 86"h; cabinet with no toe kick is 80"h

Interior configuration:

- 30 FlexCells per compartment, 60 total

Evolve 27.25"d cabinets

18.75"w cabinet

The 18.75"-wide, 27.25"-deep, 84"-high includes a counterbalanced roll-top door or solid-front or tempered glass door that can be hinged on the left or the right. Cabinet has a 4"-high toe kick.

Exterior dimensions:

- 18.75"w x 27.25"d x 84"h
- Cabinet with 6" toe kick is 86"h; cabinet with no toe kick is 80"h

Interior configuration:

- 30 FlexCells (hinged door)
- 27 FlexCells (roll-top door)

36"w cabinet

The 36"-wide, 84"-high cabinet is available with center column, no center column, or split center column. Door choices include a counterbalanced roll-top door or solid-front or tempered glass doors. Cabinet has a 4"-high toe kick.

Exterior dimensions:

- 36"w x 27.25"d x 84"h
- Cabinet with 6" toe kick is 86"h; cabinet with no toe kick is 80"h

Interior configuration:

- 30 FlexCells, 60 in cabinet with center column (hinged door)
- 27 FlexCells, 54 in cabinet with center column (roll-top door)

Product number system for Evolve FlexCell cabinets



Product Shown:
Cabinet with FlexCell, no center column
SEA2736TGNC

S	E	A	2	7	3	6	T	G			N	C
Product Line	Material	Depth	Width	Cabinet Type	Door Type	Hinge Location	Cabinet Configuration					

Product Line
SE – Evolve
Material
A – AireCore
B – AireCore with Brushed Aluminum
C – Compression Board
S – Stainless Steel
Depth (rounded to 0.25")
13 – 13"
19 – 19"
27 – 27"
Width (rounded to 0.25")
19 – 19"
27 – 27"
36 – 36"
40 – 40"

Cabinet Type
T – Tall
B – Base
U – Upper
Door Type
G – Glass
RT – Roll-top
S – Solid-front
Hinge Location
R – Right
L – Left

Cabinet Configuration
NC – No Center Column
C – Center Column
LC – Left Offset Column
RC – Right Offset Column
SC – Split Center Column
DS – Divided Shelves
AP – Accessory Pack
CC – Catheter
BC – Boxed Catheter

Products shown: Cabinet with FlexCell, no center column (SEA2736TGNC)



InnerSpace

PORTFOLIO OF PRODUCTS

Solaire Medical specializes in storage solutions for all areas of your healthcare facilities. We offer the InnerSpace® portfolio of products, a cohesive system of carts, cabinets, and modular casework that improves staff and facility efficiency and an organization's bottom line. Our InnerSpace portfolio is supplemented with the Quick line of open storage products.

When you purchase InnerSpace by Solaire Medical products, you're buying products made in the U.S.A. and from a direct manufacturer that delivers on-time and error-free. The portfolio's five-year

warranty is our promise to you that our products and service are of the highest quality.

The Quick family of open storage products have a one-year warranty.

Contact our service team with your requests or to connect with a sales representative. We're here to help you specify the best storage products for your needs, whether a standard option or customized solution. We offer complimentary storage analysis, design services, CAD support, and logistics and installation assistance.

solaire
MEDICAL