Evolve Scope Cabinet with Roll-Top Door


Evolve modular cabinets and casework create a flexible storage system that organizes a hospital's most supply-intensive departments. A clean design and tight construction offer a refined aesthetic and support infection prevention efforts.

The 19.25"-deep scope cabinet comes in 2 widths and 2 heights. It has a counterbalanced roll-top door and vented top and bottom, padded rear wall, removable drip tray, and adjustable scope managers and cord retainers. Material choices include compression board with a melamine or specified laminate surface or nonporous AireCore with Glacier White finish, Brushed Aluminum finish, or specified laminate surface. Cabinet has a 4"-high toe kick.





Evolve Scope Cabinet with Roll-Top Door

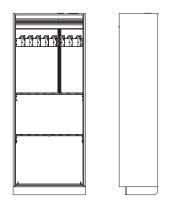

Dimensions (Dimensions are nominal)

Exterior (interior); scope hanging length:

26.75"w x 19.25"d x 60"h (25.25"w x 17"d x 47.5"h); 44" maximum 26.75"w x 19.25"d x 92"h (25.25"w x 17"d x 79.5"h); 76" maximum 36"w x 19.25"d x 92"h (34.5"w x 17"d x 79.5"h) (shown); 76" maximum 60"h cabinet with 6"h toe kick is 62"h, 56"h with no toe kick; 92"h cabinet with 6"h toe kick is 94"h, 88"h with no toe kick


Interior Configuration:

60"h cabinet holds 5 scopes 92"h x 26.75"w cabinet holds 5 or 10 scopes 92"h x 36"w cabinet holds 8 or 16 scopes

Options

- Optional Fan Ventilation Unit available
- Keyed or keyless lock options available
- Optional sloped top available
- Optional 6" toe kick or no toe kick available

Product Number Reference

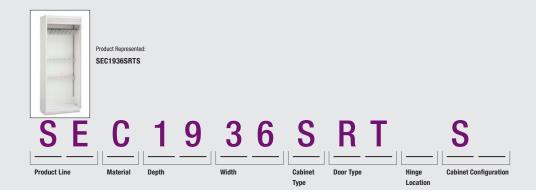

Product Numbers

26.75"w x 19.25"d x 60"h

Compression board: SEC1927SRTB AireCore: SEA1927SRTB (shown)

AireCore and Brushed Aluminum: SEB1927SRTB



26.75"w x 19.25"d x 92"h

Compression board: SEC1927SRTS; with double scope holder: SEC1927SRTD

AireCore: SEA1927SRTS (shown);

with double scope holder: SEA1927SRTD

AireCore and Brushed Aluminum: SEB1927SRTS; with double scope holder: SEB1927SRTD



36"w x 19.25"d x 92"h

Compression board: SEC1936SRTS; with double scope holder: SEC1936SRTD

AireCore: SEA1936SRTS (shown);

with double scope holder: SEA1936SRTD

AireCore and Brushed Aluminum: SEB1936SRTS; with double scope holder: SEB1936SRTD

Materials and Finishes

Cabinet shell:

Compression board: Thermally fused laminate over 0.75"-thick compression board; melamine surface or specified Wilsonart, Formica, or Pionite standard laminate

AireCore: Aluminum composite with Glacier White finish or specified Wilsonart, Formica, or Pionite standard laminate

AireCore and Brushed Aluminum finish: Aluminum composite with Brushed Aluminum finish Edge banding:

Compression board and AireCore: White ABS

AireCore with Brushed Aluminum finish: Brushed Aluminum acrylic

FlexCell panel:

Injection-molded ABS

Roll-top door and recessed channel:

Polypropylene

