



Evolve Wall-Mounted Scope Cabinet with Roll-Top Door

Evolve modular cabinets and casework create a flexible storage system that organizes a hospital's most supply-intensive departments. A clean design and tight construction offer a refined aesthetic and support infection prevention efforts.

The 26.75"-wide, 19.25"-deep, 48"-high wall-mounted scope cabinet has a counterbalanced roll-top door. It is vented on the top and bottom, has a padded rear wall, includes a removable drip tray, and has adjustable scope managers and a cord retainer. Material choices include compression board with a melamine or specified laminate surface or nonporous AireCore with Glacier White finish, Brushed Aluminum finish, or specified laminate surface.

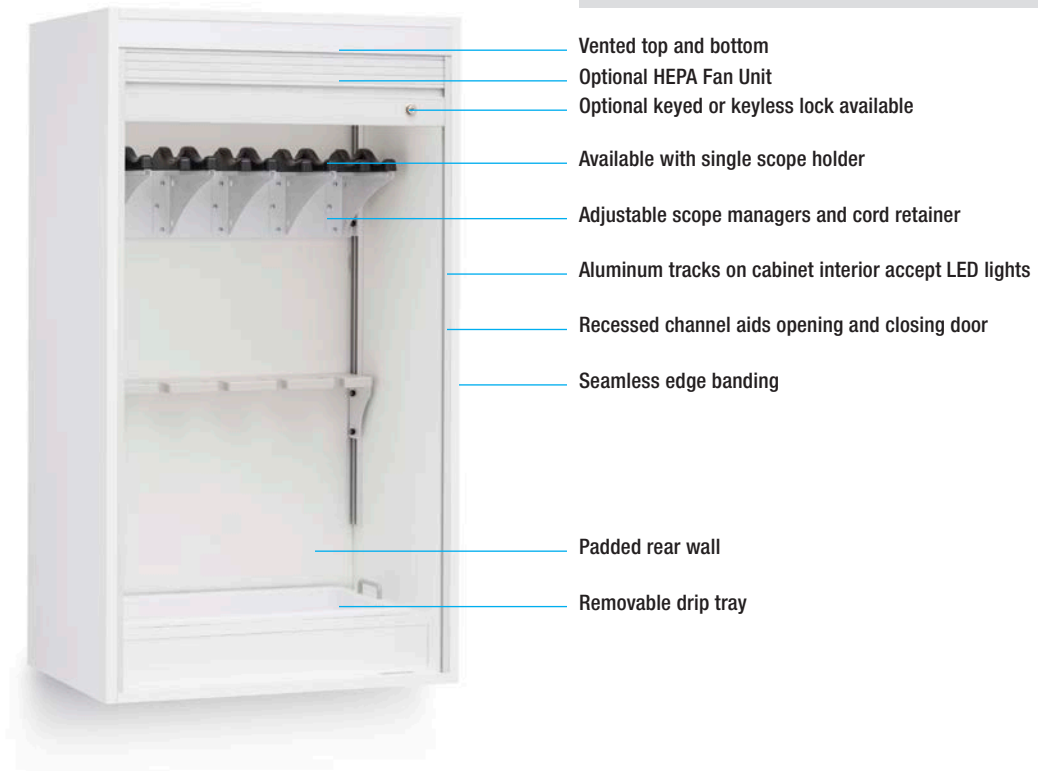

Evolve Wall-Mounted Scope Cabinet with Roll-Top Door

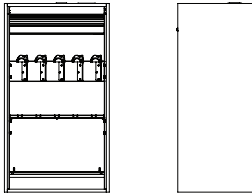

Dimensions (Dimensions are nominal)

Exterior (interior); scope hanging length:

26.75"w x 19.25"d x 48"h (25.25"w x 17"d x 39.5"h) (shown); 33" maximum

Interior Configuration:

Cabinet holds 5 scopes

Options

- Optional Fan Ventilation Unit available
- Keyed or keyless lock options available
- Optional sloped top available

Product Number Reference



Product Represented:
SEC1927SRTW

S E C 1 9 2 7 S R T W

Product Line	Material	Depth	Width	Cabinet Type	Door Type	Hinge Location	Cabinet Configuration
--------------	----------	-------	-------	--------------	-----------	----------------	-----------------------

Product Numbers



26.75"w x 19.25"d x 48"h

Compression board: SEC1927SRTW

AireCore: SEA1927SRTW (shown)

AireCore and Brushed Aluminum: SEB1927SRTW

Evolve Wall-Mounted Scope Cabinet with Roll-Top Door

Materials and Finishes

Cabinet shell:

Compression board: Thermally fused laminate over 0.75"-thick compression board; melamine surface or specified Wilsonart, Formica, or Pionite standard laminate

AireCore: Aluminum composite with Glacier White finish or specified Wilsonart, Formica, or Pionite standard laminate

AireCore and Brushed Aluminum finish: Aluminum composite with Brushed Aluminum finish

Edge banding:

Compression board and AireCore: White ABS

AireCore with Brushed Aluminum finish: Brushed Aluminum acrylic

FlexCell panel:

Injection-molded ABS

Roll-top door and recessed channel:

Polypropylene



Storage solutions for healthcare

Solairé Medical | P.O. Box 2808 | Grand Rapids, Michigan 49501
888-435-2256 phone | 616-743-5944 fax | solairemedical.com

© 2018 Solairé Medical, Grand Rapids, MI 49501
© InnerSpace and Scope Manager are registered trademarks of Solairé Medical, Grand Rapids, MI 49501
0918